


**RUSCH**

**CLOSED URINARY  
DRAINAGE SYSTEMS**

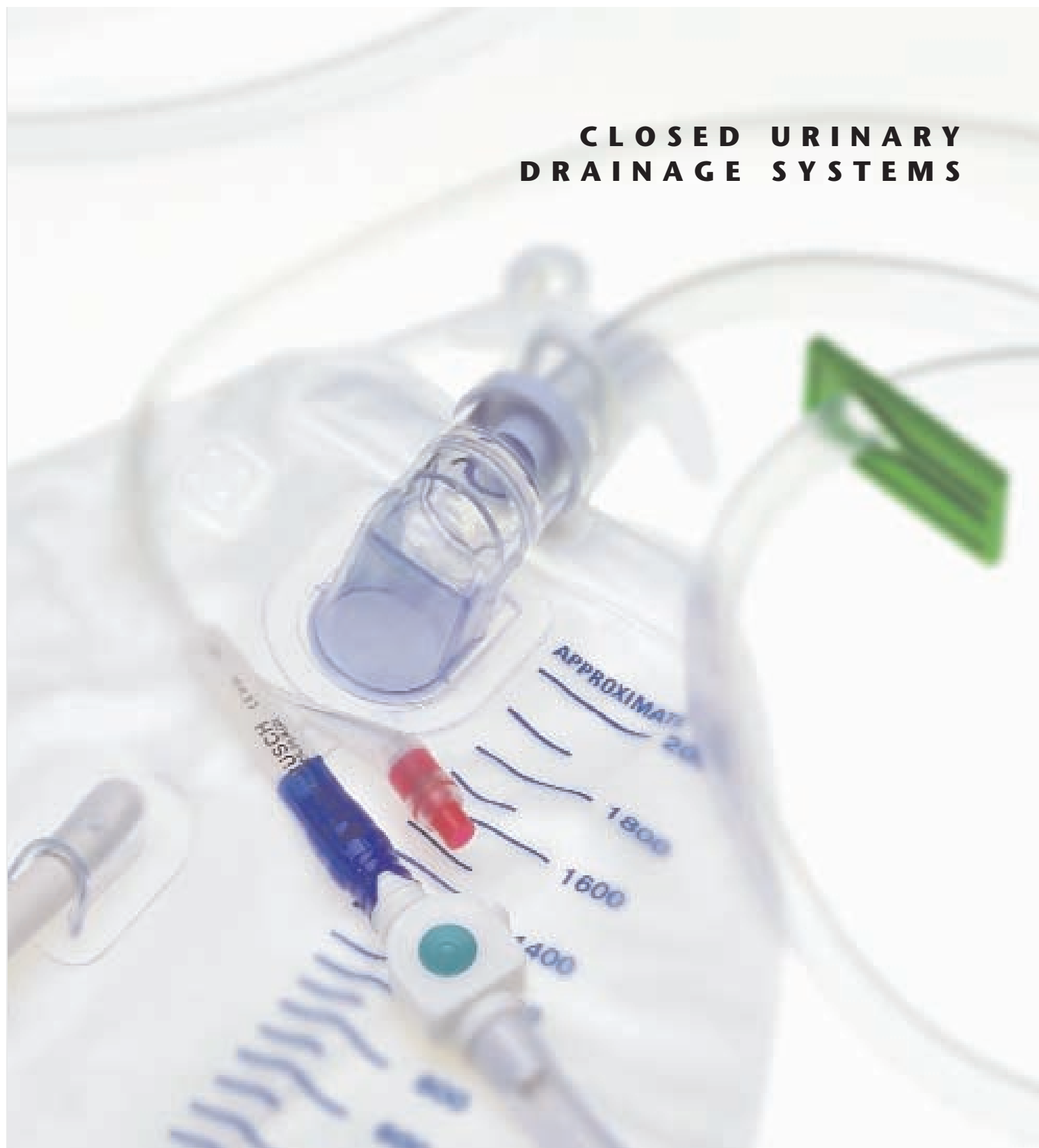

FOR MAXIMUM  
PROTECTION OF THE  
PATIENT AGAINST  
NOSOCOMIAL  
URINARY TRACT  
INFECTIONS

## FOR MAXIMUM PROTECTION OF THE PATIENT AGAINST NOSOCOMIAL URINARY TRACT INFECTIONS

About 90% of all nosocomial urinary tract infections are caused by urinary drainage. Apart from the health risks for the patient, such infections also present a big economic problem.

Rüsch offers a variety of products which reduce the risk of infection. The closed urinary drainage system represents one of the most important elements of prophylaxis against infections.


Rüsch urinary drainage systems meet all the hygienic demands which are required for a high-quality closed system for long-term medical care.

### FLEXIBLE AND RESISTANT TO KINKING


Large-lumen transparent  
tube with sliding clamp  
and sheet clip.

### HYGIENIC AND QUICK EMPTYING

Non-dripping bottom outlet  
with one-hand clamp.  
Doghouse holder to fix the outlet  
tube, to prevent contamination,  
and to achieve additional floor  
clearance.

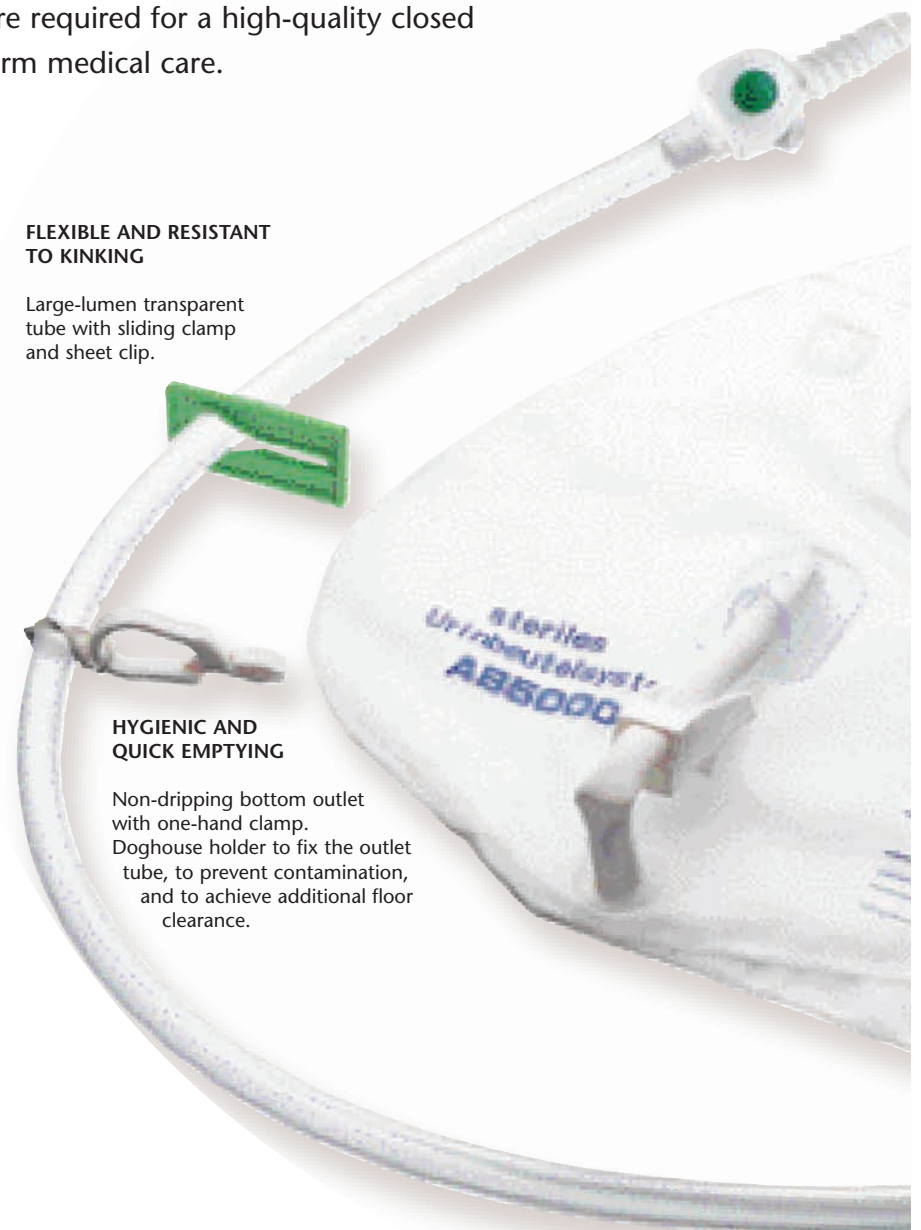

### URINE SAMPLING – SAFE, EASY AND WITHOUT THE NEED FOR A NEEDLE

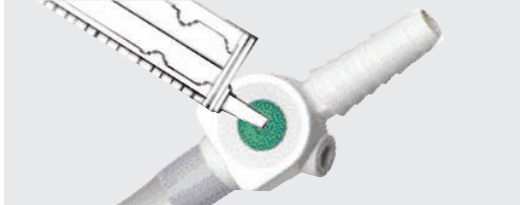


Urine sampling is carried out at the ridged universal connector with a Luer syringe without a needle.

Thus, needle stick injuries are eliminated and the risk of invading bacteria is greatly reduced.

### SAFE AND SOLID CONNECTION

Vented universal ridged connector with tapered cone-like steps – suitable for all common catheters. AB 5000 also preconnected with silicone balloon catheter.

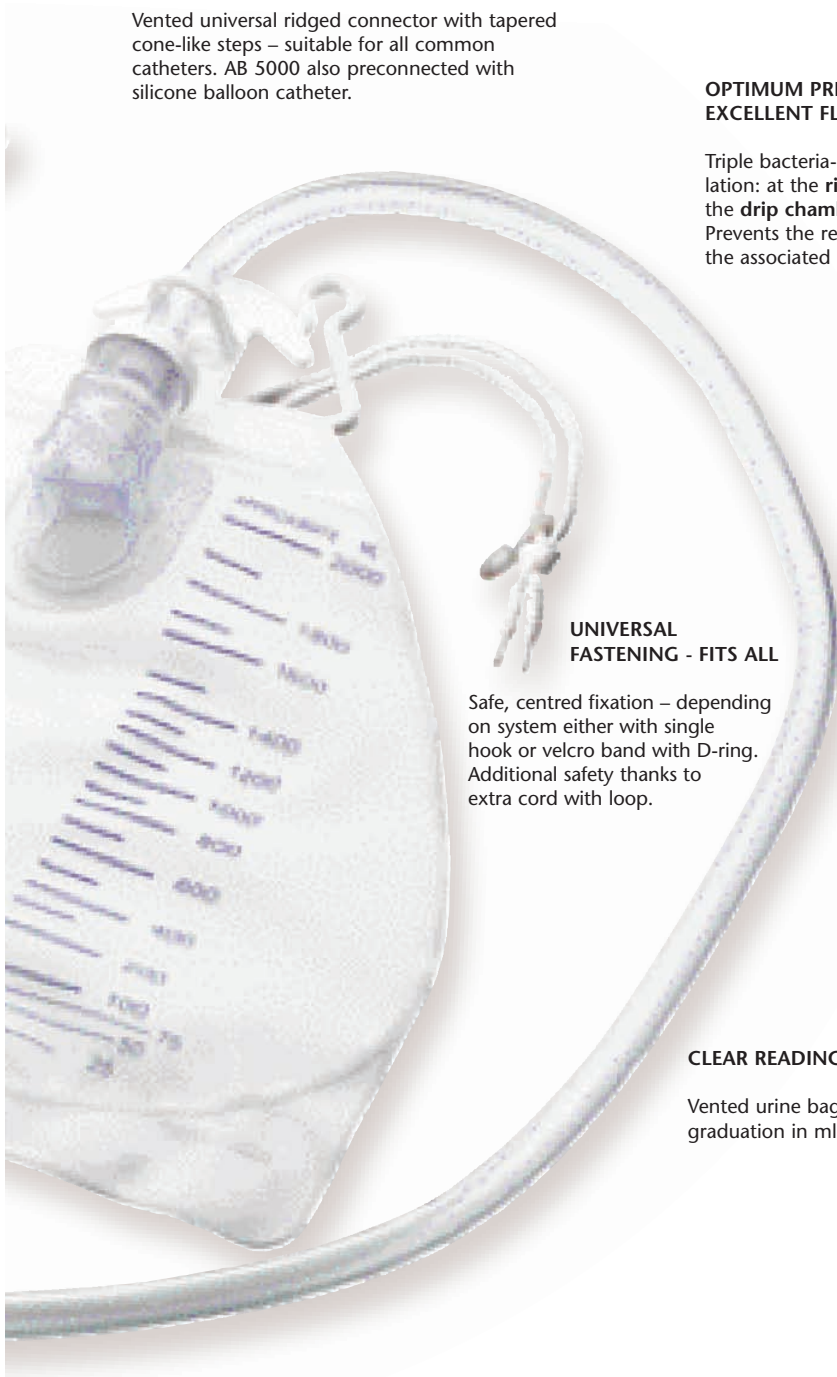

### OPTIMUM PRESSURE BALANCE AND EXCELLENT FLOW PERFORMANCE

Triple bacteria-proof and liquid-tight ventilation: at the **ridged universal connector**, the **drip chamber** and the **urine bag**. Prevents the reflux of urine in the tube and the associated risk of infection.

### UNIVERSAL FASTENING - FITS ALL

Safe, centred fixation – depending on system either with single hook or velcro band with D-ring. Additional safety thanks to extra cord with loop.

### CLEAR READING


Vented urine bag with clear graduation in ml.

### VENTED DRIP CHAMBER WITH INTEGRATED CHECK VALVE

*The vented drip chamber with high drip level and safely functioning diaphragm valve reliably prevents retrograde rising infections and urine reflux. Optimal pressure compensation and safe outlet performance thanks to bacteria-proof ventilation. Safe even under high pressure load.*



## PRECONNECTED

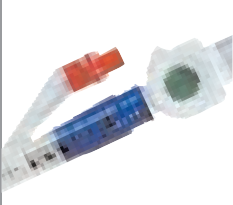

**PRECONNECTED:**  
MAXIMUM PROPHYLAXIS AGAINST INFECTIONS DURING STANDARD MEDICAL CARE

*Reduction of intraluminally rising urinary tract infections: The connection is safeguarded with a protective foil – thus preventing unintentional disconnection with the associated invasion of bacteria. (If necessary, the foil can be removed)*

*Economic and quick: No need for connection – saves time and reduces the number of application mistakes during catheterisation.*

*The pre-printed adhesive labels also save time during documentation – manual notes are made redundant.*

### Closed urinary drainage system with silicone balloon catheter for standard drainage under long-term conditions

The positive characteristics of the AB 5000 and the Rüsç Brilliant silicone balloon catheter combined with the hygienic advantages of a preconnected system: maximum safety, high patient comfort and economic system solution.



#### Preconnected

Sterile one-way product, latex-free, Qty 20 consists of: **Ref. 17 11 00**

##### Rüsç Brilliant silicone balloon catheter

cylindrical closed tip **Size (Ch.)**  
2 opposed eyes **12 / 14 / 16 / 18**  
length approx. 41 cm, balloon 10 ml

##### Closed urinary drainage system AB 5000

capacity 2000 ml  
length of tube 120 cm

**SilFlate - prefilled syringe** **Ref. 85 05 50\***

capacity 10 ml

#### AB 5000:

- **Urine sample port without needle** at the ridged universal connector
- **Drip chamber** with integrated check valve
- **Triple bacteria-proof and liquid-tight ventilation** at the ridged universal connector, drip chamber and urine bag
- **Clear graduation**
- **Flexible and non-kinking tube** with sliding clamp and sheet clip
- **Optimum fastening of the system** by means of a single hook and an additional cord with loop
- **Non-dripping bottom outlet** for one-hand operation, with doghouse holder

#### Rüsç Brilliant silicone:

- **Transparent 2-way catheter** with radiopaque tip
- **Silicone** guarantees constant elasticity and biocompatibility
- **Large eyes and large inner lumen** permit high drainage performance
- **Minimal tendency to incrustation** thanks to smooth surfaces
- **Atraumatic tip and eyes**
- **Enclosed prefilled syringe** for inflation of the balloon, in sterile packaging, 10 ml sterile 10% aqueous glycerine solution



\* to be ordered separately

## AB 5000

- **Urine sample port without needle** at the ridged universal connector
- **Drip chamber** with integrated check valve
- **Triple bacteria-proof and liquid-tight ventilation** at the ridged universal connector, drip chamber and urine bag
- **Clear graduation**
- **Flexible and non-kinking tube** with sliding clamp and sheet clip
- **Optimum fastening of the system** by means of a single hook and an additional cord with loop
- **Non-dripping bottom outlet** for one-hand operation, with doghouse holder

Closed urinary drainage system for standard drainage under long-term conditions



### AB 5000

Capacity	Length of tube	Qty	Ref.
2000 ml	120 cm	20	67 30 41

Sterile one-way product, latex-free

## AB 6000

- **Mounted double measurement chamber** with integrated measurement function for small amounts
- **Urine sample port without needle** at the ridged universal connector
- **Drip chamber** with integrated check valve
- **Triple bacteria-proof and liquid-tight ventilation** at the ridged universal connector, drip chamber and urine bag
- **Clear graduation**
- **Flexible and non-kinking tube** with sliding clamp and sheet clip
- **Optimum fastening of the system** by means of a velcro band with D-ring and an additional cord with loop
- **Non-dripping bottom outlet** for one-hand operation, with dog house holder

Closed urinary drainage system with urine meter (double measurement chamber), to check the hourly discharge of urine, also under long-term use



### AB 6000

Capacity	Length of tube	Qty	Ref.
2000 ml	150 cm	10	67 30 61

Sterile one-way product, latex-free

AB 6000: MOUNTED  
PIVOTING DOUBLE  
MEAS-



UREMENT CHAMBER

1ml-precise graduation up to 70 ml for hourly measurement (with sampling port for fresh urine). Overflow into a 450 ml measurement chamber for larger intervals. Easy and quick emptying of the measurement chamber without the removal of any parts.

## RÜSCH TUR BAG



Closed urinary drainage system for post-operative bladder irrigation

### FIRM FIXATION AT BEDSIDE

Safe, centred positioning with single hook – fits all common bedrails. For extra safety equipped with an additional cord with loop.

Exact reading thanks to balanced position of bag.

- **Urine sample port without needle** at the ridged universal connector
- **Drip chamber** with integrated check valve
- **Triple bacteria-proof and liquid-tight ventilation** at the ridged universal connector, drip chamber and urine bag
- **Clear graduation**
- **Flexible and non-kinking tube** with sliding clamp and sheet clip
- **Optimum fastening of the system** by means of a single hook and an additional cord with loop
- **Non-dripping bottom outlet** for one-hand operation, with doghouse holder

The Rüsç TUR Bag for postoperative care. The large capacity of 4000 ml considerably reduces the emptying frequency, especially with irrigation after surgical interventions.

### Rüsç TUR Bag

Capacity	Length of tube	Qty	Ref.
4000 ml	120 cm	20	67 30 51

Sterile one-way product, latex-free


## BELLY BAG®




**Belly Bag® for all mobile patients who are discharged with a bladder catheter**

The first urine bag which can be worn like a bum bag – wearer comfort day and night.

- **Soft fabric belt** ensures that the bag remains securely in place under the patient's clothes
- **Skin-friendly and soft reverse side** ensures high level of wearer comfort
- **Non-dripping bottom outlet with twist-off tap** enables hygienic and quick emptying of the bag
- **Handy retaining band** accommodates the bottom outlet
- **Ridged universal connector** for all common catheters
- **Check valve:** integrated anti-reflux valve at the ridged universal connector
- **Urine sample port** at the ridged universal connector (only Ref. B 1000 P)
- **Clear graduation**

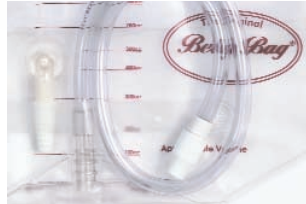

Model B 1000 CT

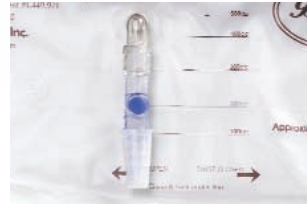

Model B 1000 P

Function: Due to the pressure\* of the bladder muscles, urine is reliably transported via the catheter out of the bladder into the bag.

\* The residual pressure of the bladder amounts to 10 – 25 cm H<sub>2</sub>O. Only 6 cm H<sub>2</sub>O are required for urine transport from the bladder into the bag.

### Belly Bag®

Description	Ref.
1000 ml	<b>B 1000</b>
1000 ml, with urine sample port	<b>B 1000 P</b>
1000 ml, with approx. 60 cm tube	<b>B 1000 CT</b>
Spare belt	<b>PO 337</b>
Sterile one-way product, latex-free, Qty 10	



**BELLY BAG®:**  
COMFORTABLE SOLUTION FOR ALL MOBILE PATIENTS

Considerable improvement of mobility and thus of the patient's quality of life – discreet and simple to hide under your clothes. No change from leg to bedside bag necessary any more.

Reduced risk of unintentional catheter extraction, as there is no drainage tube.

Necrosis of the urethra, as caused by hanging urine bags with men, are avoided with the Belly Bag®.

---

**RÜSCH**

---

 **HUDSON RCI**

---

**RÜSCH CARE**

---

**MEDICAL**  **SERVICE**

---

**WECK**

---

**Pilling**

---

**PLEUR-EVAC**

---



---

**TAUT**

---

Teleflex Incorporated is a diversified company with business interests in the aerospace, automotive, marine, industrial, and medical markets. Teleflex Medical has been for years a global supplier of innovative and quality products and services. Teleflex Medical is comprised of three divisions or business units: Surgical, Medical and OEM.

[www.ruesch.de](http://www.ruesch.de)  
[www.teleflex.com](http://www.teleflex.com)

**GERMANY** TELEFLEX MEDICAL GmbH, P.O. Box 1180, 71385 Kernen  
Phone +49 (0)7151 406 0, Fax +49 (0)7151 406 150, [info.de@teleflexmedical.com](mailto:info.de@teleflexmedical.com)

**AUSTRIA** RÜSCH AUSTRIA Gesellschaft m.b.H., A-1090 Wien  
Phone +43 (0)1 402 47 72, Fax +43 (0)1 402 47 72 77, [rueschaustria@teleflexmedical.com](mailto:rueschaustria@teleflexmedical.com)

**BELGIUM** TELEFLEX MEDICAL BE sprl, B-1070 Brussels, Phone +32 (0)2 333 24 60  
Fax +32 (0)2 332 27 40, [info.be@teleflexmedical.com](mailto:info.be@teleflexmedical.com)

**CANADA** RUSCH CANADA, Oakville - Ontario L6M 3E3  
Phone +1 (0)905 469 8489, Fax +1 (0)905 469 8490

**FRANCE** TELEFLEX MEDICAL SAS, F-31460 Le Faget, Phone +33 (0)5 62 18 79 40  
Fax +33 (0)5 61 83 35 84, [info.teleflexmedical.france@teleflexmedical.com](mailto:info.teleflexmedical.france@teleflexmedical.com)

**INDIA** TELEFLEX MEDICAL Pvt. Ltd., Chetpet, Chennai - 600 031  
Phone +91 (0)44 2836 5040/-41, Fax +91 (0)44 2836 0682, [mktg@teleflex-ind.com](mailto:mktg@teleflex-ind.com)

**ITALY** TELEFLEX MEDICAL S. r. l., 20039 Varedo (MI), Phone +39 0362 58911  
Fax +39 0362 5891 888, [info.it@teleflexmedical.com](mailto:info.it@teleflexmedical.com)

**JAPAN** PILLING WECK (ASIA) PTE. LTD., Minato-ku, Tokyo 105-0001  
Phone +81 (0)3 5404 3415, Fax +81 (0)3 5404 3401

**NETHERLANDS** TELEFLEX MEDICAL B.V., P.O. Box 40, 3640 AA Mijdrecht,  
Phone +31 (0)297 540 703, Fax +31 (0)297 533 726 / 540 166, [info.nl@teleflexmedical.com](mailto:info.nl@teleflexmedical.com)

**SINGAPORE** PILLING WECK (ASIA) PTE. LTD., Lavender Place  
Singapore 338750, Phone +65 (0)6294 1019, Fax +65 (0)6299 6132

**SPAIN** TELEFLEX MEDICAL S.A., 03530-La Nucía (Alicante)  
Phone +34 (0)96 687 51 01/-02/-03, Fax +34 (0)96 587 35 69, [info.es@teleflexmedical.com](mailto:info.es@teleflexmedical.com)

**SWITZERLAND** RÜSCH (SCHWEIZ) AG, CH-3123 Belp  
Phone +41 (0)31 819 73 90, Fax +41 (0)31 819 73 87, [info@ruesch-schweiz.ch](mailto:info@ruesch-schweiz.ch)

**UNITED KINGDOM** TELEFLEX MEDICAL, High Wycombe, Bucks. HP12 3ST  
Phone +44 (0)1494 53 27 61, Fax +44 (0)1494 52 46 50, [info.uk@teleflexmedical.com](mailto:info.uk@teleflexmedical.com)

**USA** TELEFLEX MEDICAL, Durham, NC 27709  
Phone +1 (800) 234 9325 , Fax +1 (800) 932 5329, [csrusch@teleflexmedical.com](mailto:csrusch@teleflexmedical.com)