

AUTOCLAVABLE ELECTRODES

All the PRECISION You Require in Surgery with the VERSATILITY You Need.

All Ellman® electrodes are designed for use in high frequency (4 MHz) RF Generators. The insulated shafts can be bent or contoured to facilitate use in various surgical procedures.



Cat. No. XA



Cat. No. XB



Cat. No. XC



Cat. No. XD



Cat. No. XT



Cat. No. XP


CHOICE OF WIRE GAUGE
 — Fine - Used when minimal tissue alteration is desired
 — Regular - Used to produce more coagulation

	<u>Cat. No.</u>
	TA1
	TA2
	TA9
	TA3
	TF1
	TA8
	TBB
	TB1, TB4
	TB9, TB10
	TB2, TB5
	TB3, TB6
	TD3
	TD8
	TD9
	TD10
	TD11
	TC1, TC4
	TC3, TC6
	TC17, TC18
	TC19, TC20
	TC21, TC22
	TC9, TC10
	TC7, TC8
	TC11, TC12
	TC13, TC14
	TC15, TC16
	TP1, TP2
	TP3, TP4
	TP5, TP6
	TP7, TP8
	TP9, TP10

AUTOCLAVABLE ELECTRODES

Needle Electrodes (A Series)

- Ideal for incisions & excisions in precise surgical procedures
- Choice of gauge is dependent on degree of coagulation desired

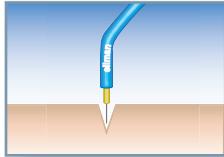


Round Loop Electrodes (B Series)

- Ideal for shaving lesions and contouring tissue
- Excellent for biopsy




Diamond Loop Electrodes (C Series)

- Provides one-stroke elliptical excision

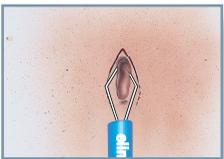

Ball Electrodes (D Series)

- Ideal for coagulation
- Heat seals blood vessels rather than burning tissue
- Minimal adjacent tissue destruction

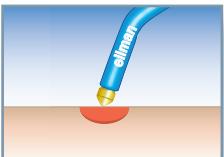

Broad Needle Electrodes

- Ideal for Fulguration and Desiccation procedures

Oval Loop Electrodes (P Series)

- Ideal for shaving lesions and contouring tissue
- Alternative to a round loop when working in tight spaces


Triangle Loop Electrodes (T Series)

- Ideal for shaving lesions and contouring tissue

Vari-Tip™ Microdissection Needle


PRECISION:

Vari-Tip wire can be shortened, resulting in a superficial incision, by pushing against your thumbnail. (i.e. skin incisions.)

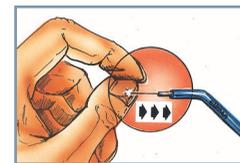

VALUE:

Vari-Tip wire may be extended and cut for continued use. When the wire has been completely trimmed down, simply insert the extra wire provided.

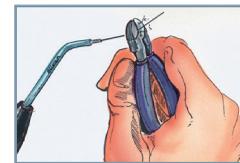

Electrode Block

ADVANTAGES

- Organizes electrode tips for simple identification and selection
- Holds 18 electrodes upright
- Durable lucite




Cat. No.	Description
H1	Electrode Block (with 18 holes) 1/16"
QH1	3/32" (18 holes)
H1A	Combination with both shaft sizes (18 holes)

Retractable Electrode Protector Stand

ADVANTAGES

- Protects electrodes from damage and breakage
- Presents electrodes for easy identification and access
- Holds both 1/16" and 3/32" shafts
- Holds 12 electrodes

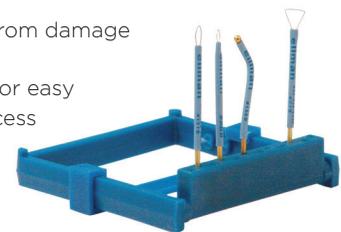

Cat. No.	Description
H66	Small (2" Tall)
H66A	Large (4" Tall)

See our Accessories Catalog for our line of specialty electrodes and disposable electrodes. Contact our technical sales representatives for shipping costs, delivery options, and discounted volume packages.

Ellman • 400 Karin Lane • Hicksville, NY 11801 • Tel (516) 594.3333 • www.ellman.com

Cynosure and Ellman are registered trademarks and Vari-Tip is a trademark of Cynosure, Inc.

© 2015 Cynosure, Inc • CC11131 REV B